

# Urostomy pouches

Clearly the natural choice.



**OakMed** 

Giving you more choice.

# All Oakmed urostomy pouches benefit from the following features:

---

Available in four wafer types:

- **Hydrocolloid**
- **Microskin®**
- **Bi Form®**
- **Soft Convex**



## Unique figure-of-eight shape

Provides maximum capacity with minimum bulk.


## Non-return valve

All Oakmed urostomy pouches have a non-return valve inside to stop the urine back-flowing over the stoma and causing infections. Stents can be easily inserted through the non-return valve or curled into the top section of the pouch.




## Sizes

All four flange types are available in a large range of pre-cut sizes to suit the majority of urostomies.

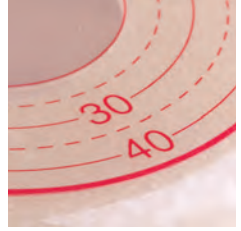

## Soft fabric cover

Quiet, non-cling fabric for maximum comfort.

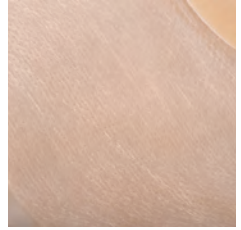

## Tap

Easy to empty bung, safe from leaks when closed.




## Universal connector

Convenient and easy to use bung outlet, this can be used with any night drainage bag with the connector provided.



## Convenient split fabric front cover

For accurate positioning and easy observation. Hydrocolloid pouch also available without soft front cover.

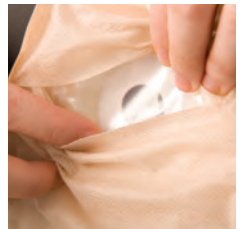



## **This Hydrocolloid is hypo-allergenic and skin friendly.**

The Oakmed Hydrocolloid has been specifically designed with an extremely flexible tapered edge which is integral to the pouch, this eliminates the need for flange extensions making pouch changes quicker and less complex.

- The flange will conform easily to all skin creases and folds giving added security and confidence.
- The Hydrocolloid is of an optimum thickness specifically for urostomists, this helps to prevent breakdown of the adhesive by the urine which can lead to reduced wear time and subsequent leakages.
- The Hydrocolloid has a 'high-tack' adhesion, designed to help prevent urine tracking behind the wafer. Again this allows for added security.
- The Hydrocolloid is soft and flexible and can be removed easily.
- Available with and without soft cover.
- Cuts up to 60mm and is available in a range of pre-cut sizes, (clear pouch pre-cut with 13mm starter hole only).



**Unique to Oakmed the Microskin® adhesive is thin, transparent and flexible, moulding easily to all body contours and creases, significantly reducing the risk of leakage.**



Microskin® offers an all new level of comfort and security.

- Barely visible on your skin when it is in place.
- Moves with you like a second skin.
- Because it is transparent, allows you to visualise the health of your skin underneath the adhesive, and also check for any leakages occurring.
- Breathable and moisture vapour permeable for healthier skin.
- Not affected or dissolved by urine allowing for a longer wear time and added security.
- Supports the pouch without the need for a belt or additional flange extensions.
- Waterproof and will stay in place whilst showering, bathing or swimming. Towel dry as normal.
- Cuts up to 50mm and is available in a range of pre-cut sizes.

## TCI300-mini Urostomy Microskin®

- This pouch has all the benefits of the Microskin® adhesive.
- Available in a starter hole of 13mm, which cuts up to 50mm.
- Ideal for urostomists to use for leisure purposes, i.e. swimming.
- Has a clear front and a soft fabric backing for added comfort.
- Easy to use tap which can be connected to a drainage system.
- Can also be used on small wound sites.



Bi Form®

biform®

**Bi Form® combines the benefits of both  
Microskin® and Hydrocolloid adhesives.**

Hydrocolloid protects the peri-stomal skin and is of an optimum thickness to help prevent breakdown by urine.

Microskin® outer gives added security and confidence and will conform to all skin creases and folds. It also eliminates the need for any flange extension products.

Microskin® is clear on the skin which allows for early detection of any leakage.

Ideal for use with para-stomal hernias as Microskin® will conform around the hernia and give an excellent seal, for added security against leakage.

Hydrocolloid around the peri-stomal skin makes this product ideal for inserting stents in the post-operative period.

Cuts up to 45mm and is available in a range of pre-cut sizes.



Bi Form®



# Hydrocolloid Soft Convex



**The Hydrocolloid Soft Convex flange is low profile whilst also being substantial enough to be effective against leakage, i.e. with flush or retracted stomas.**

The Hydrocolloid is of an optimum thickness to help prevent breakdown of the flange by the urine. This gives added confidence against leakage.


The wafer has a unique, tapered, integral border. This helps the flange to mould over creases and dips easily and eliminates the need for flange extensions.

Providing a 'high-tack' superior adhesion, this wafer gives a superb fit over the body contours.

The Hydrocolloid is gentle on the skin, non-irritant and is easily removed.

The Soft Convex pouch has convenient belt loops for those who prefer the security of a belt.

Cuts to a maximum of 44mm and is available in a range of pre-cut sizes.

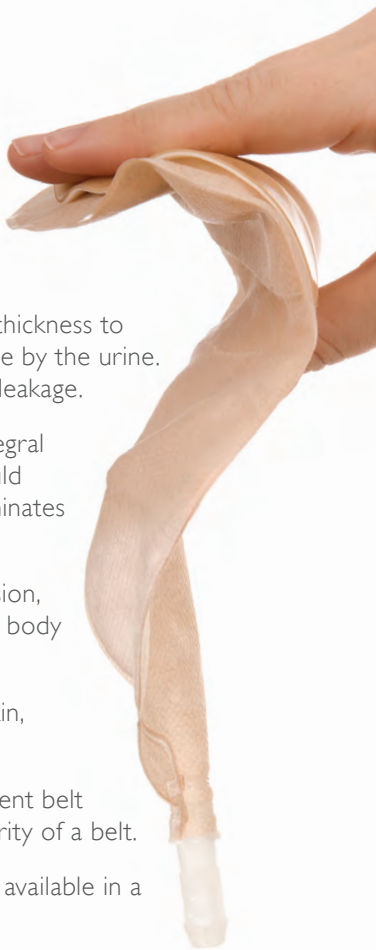

# Soft Convex



Option Hydrocolloid Soft Convex Urostomy 25mm Plateau		Packs of 20
Hole Size	Plateau size	Window
Cut to fit 15mm to 24mm	25mm	SCU25-0015
20mm	25mm	SCU25-0020
Option Hydrocolloid Soft Convex Urostomy 35mm Plateau		Packs of 20
Hole Size	Plateau size	Window
Cut to fit 15mm to 34mm	35mm	SCU35-0015
20mm	35mm	SCU35-0020
25mm	35mm	SCU35-0025
30mm	35mm	SCU35-0030
Option Hydrocolloid Soft Convex Urostomy 45mm Plateau		Packs of 20
Hole Size	Plateau size	Window
Cut to fit 15mm to 44mm	45mm	SCU45-0015
20mm	45mm	SCU45-0020
25mm	45mm	SCU45-0025
30mm	45mm	SCU45-0030
35mm	45mm	SCU45-0035
40mm	45mm	SCU45-0040

SCU25 SCU35 SCU45 Opaque Soft Cover Both Sides Window

Mini Urostomy Microskin with window & soft cover to both sides		Packs of 20
Hole Size		Microskin
Cut to fit 13mm - 50mm		TCI300

Option Urostomy with window & soft cover to both sides	Packs of 20	Packs of 20	Packs of 20
Hole Size	Bi Form	Microskin	Hydrocolloid
Starter hole 13mm See below	BF-1300k	MJ1300	GC1300
16mm	BF-1600k	MJ1600	GC1600
19mm	BF-1900k	MJ1900	GC1900
22mm	BF-2200k	MJ2200	GC2200
25mm	BF-2500k	MJ2500	GC2500
32mm	BF-3200k	MJ3200	GC3200
38mm	BF-3800k	MJ3800	GC3800

Option Urostomy Hydrocolloid clear soft cover to one side		Packs of 20
Hole Size		Hydrocolloid Clear
Starter hole 13mm See below		GC1300C
BF-1300k cuts up to 45mm		
MJ1300 cuts up to 50mm		
GC1300 cuts up to 60mm		

**OakMed**   
Giving you more choice.

Unit B, The I/O Centre, Barn Way, Lodge Farm Industrial Estate, Northampton, NN5 7UW  
Freephone: 0800 592 786 | Main: 01604 586 529 | Fax: 01604 584 784 | [www.oakmed.co.uk](http://www.oakmed.co.uk)

