



## Alginate wafer

For skin that is prone to soreness, clinically proven to stimulate healing.




## New filter system

Eliminates odour and spotting. Also helps to reduce ballooning and pancaking.



## Sizes

A comprehensive range of cut to fit and pre-cut sizes.

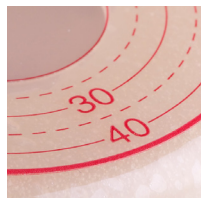

## Hygienic self-sealing closure

Convenient and washable.



## Closed option

A range of closed pouches in cut to fit and pre-cut sizes.



Alginate Soft Convex Pouches		Drainable	Closed
Option Alginate Soft Convex 25mm Plateau		Packs of 20	Packs of 20
Hole Size	Plateau size	Window	Window
Cut to fit 20mm to 24mm	25mm	SCA25-4120KV	SCA25-0420K
Option Alginate Soft Convex 35mm Plateau		Packs of 20	Packs of 20
Hole Size	Plateau size	Window	Window
Cut to fit 20mm to 34mm	35mm	SCA35-4120KV	SCA35-0420K
25mm	35mm	SCA35-4125KV	SCA35-0425K
30mm	35mm	SCA35-4130KV	SCA35-0430K
Option Alginate Soft Convex 45mm Plateau		Packs of 20	Packs of 20
Hole Size	Plateau size	Window	Window
Cut to fit 20mm to 44mm	45mm	SCA45-4120KV	SCA45-0420K
25mm	45mm	SCA45-4125KV	SCA45-0425K
30mm	45mm	SCA45-4130KV	SCA45-0430K
35mm	45mm	SCA45-4135KV	SCA45-0435K
40mm	45mm	SCA45-4140KV	SCA45-0440K

SCA25 SCA35 SCA45 Opaque, Filter, Soft Cover to Both Sides, Window

Hydrocolloid Soft Convex Pouches		Drainable	Closed
Option Hydrocolloid Soft Convex 25mm Plateau		Packs of 20	Packs of 20
Hole Size	Plateau size	Window	Window
Cut to fit 20mm to 24mm	25mm	SC25-4120KV	SC25-0420K
Option Hydrocolloid Soft Convex 35mm Plateau		Packs of 20	Packs of 20
Hole Size	Plateau size	Window	Window
Cut to fit 20mm to 34mm	35mm	SC35-4120KV	SC35-0420K
25mm	35mm	SC35-4125KV	SC35-0425K
30mm	35mm	SC35-4130KV	SC35-0430K
Option Hydrocolloid Soft Convex 45mm Plateau		Packs of 20	Packs of 20
Hole Size	Plateau size	Window	Window
Cut to fit 20mm to 44mm	45mm	SC45-4120KV	SC45-0420K
25mm	45mm	SC45-4125KV	SC45-0425K
30mm	45mm	SC45-4130KV	SC45-0430K
35mm	45mm	SC45-4135KV	SC45-0435K
40mm	45mm	SC45-4140KV	SC45-0440K

SCA25 SCA35 SCA45 Opaque, Filter, Soft Cover to Both Sides, Window



Giving you more choice.

Unit B, The I O Centre, Barn Way, Lodge Farm Industrial Estate, Northampton, NN5 7UW  
 Freephone: 0800 592 786 | Main: 01604 586 529 | Fax: 01604 584 784 | [www.oakmed.co.uk](http://www.oakmed.co.uk)



# Soft convex pouches

Naturally, the flexible choice.



**OakMed** 

Giving you more choice.



### Hydrocolloid wafer

Soft and flexible with a unique tapered border.



### Unique low profile

Effective against leakage.




### Soft fabric cover

Quiet, non-cling fabric for maximum comfort.


## Convenient split fabric front cover

For accurate positioning and easy observation.

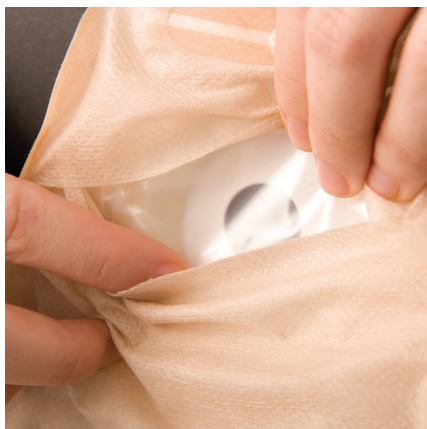

### Extended outlet

For ease of emptying.

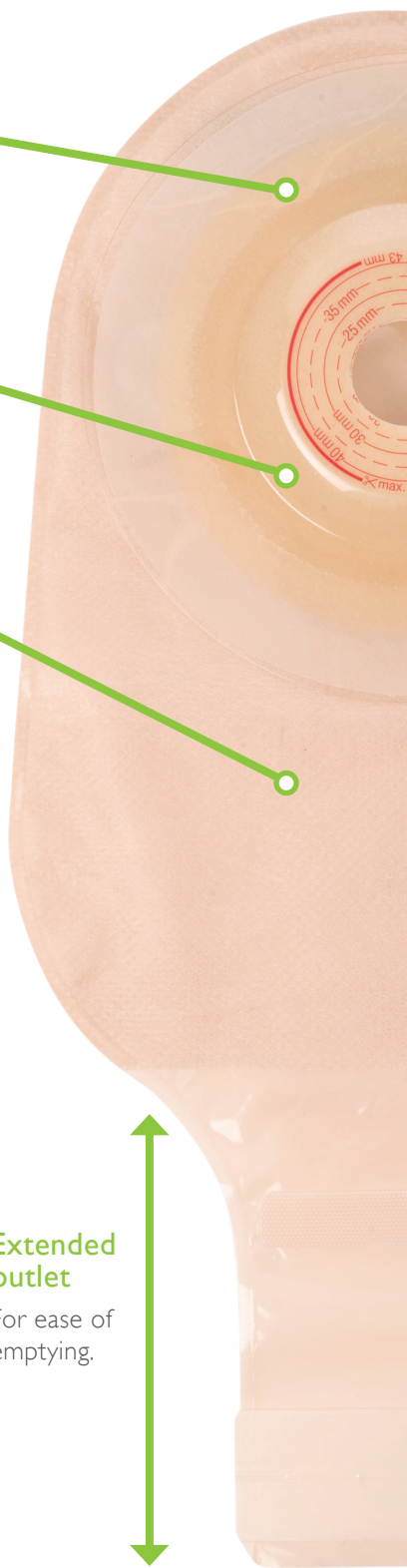

## Maximum confidence

---

### **The advantages of our new soft convex wafers are many:**

- As you move it moves with you.
- It won't crease or lift on the edges, meaning less leakage and the pouch stays on longer.
- Fits the contours of your body, so more comfortable to wear.
- More confidence to enjoy a normal life.

## Unique low profile shape

---



The convex flanges are low profile whilst also being substantial enough to be effective against leakage with flush or retracted stomas.

**Consult a Stoma Care Specialist before using Convex pouches.**

Flexible





**A choice of two superior skin protectors.**

## Alginate

The Alginate wafer has both preventative and curative properties. Alginates are natural derivatives of brown seaweed and are used extensively in wound care and have two main properties:



1. They are clinically proven to initiate skin healing and prevent tissue breakdown.
2. They also absorb any moisture, ie. bleeding or serous fluid whilst maintaining the healing process.


## Hydrocolloid

Incredibly soft and flexible, this Hydrocolloid wafer has a unique tapered border which is integral to the wafer.

Providing superior adhesion, this wafer gives a superb fit over the body contours.

**Both adhesives are slim and easy to bend when fitting the pouch over the stoma.**

### MPC® Filter - Moisture Protected Composition

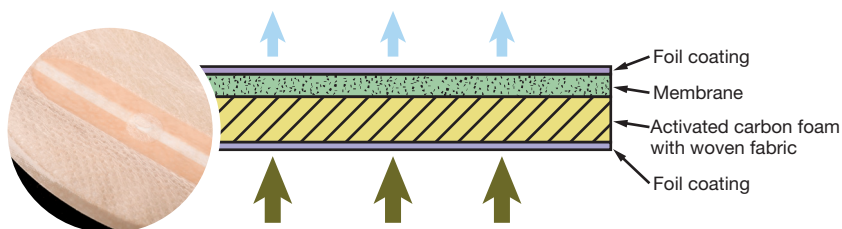

## Filter

Completely new filter system with a Moisture Protected Composition (MPC) for greater confidence and comfort. The structure of the filter allows gas to pass through, but not liquid. This eliminates odour and any spotting through the filter onto clothing. It also helps to reduce ballooning and pancaking.



## Soft fabric cover

The cover is made of soft, quiet fabric that does not cling to clothing, giving maximum comfort.

## Easy to use drainable pouch




The drainable pouch has been designed with an extended outlet which makes emptying and positioning easier. The closure tab is situated higher up the outlet for more effective cleaning and can be used to aid opening of the outlet. The tab is washable and can be folded under the soft cover for added discretion.



Comfortable

## Additional choices in both Alginate and Hydrocolloid.

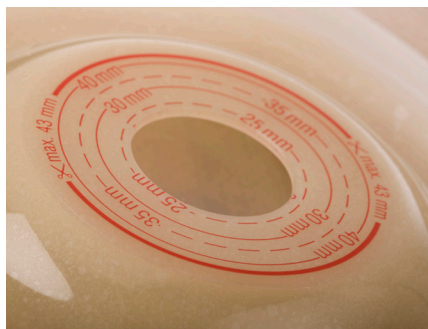

### Sizes

---

The convex range of pouches are available in 3 plateau sizes along with a range of pre-cut sizes to suit the majority of stomas.

### Soft Convex Closed

---

We offer a comprehensive range of closed pouches using the same wafers and with the same features and benefits as the drainable pouches.



### Soft Convex Drainable

---

The Soft Convex Drainable is a clipless pouch that provides maximum security and confidence. With its unique features it is comfortable and easy to use.

### Peace of mind

---

The convex range pouches all have soft flexible discreet belt loops for added security, for those who prefer the option of a belt.

The outlet on the drainable pouch is suitable for use with a separate clip if preferred.

The Alginate adhesive with its natural skin healing properties eliminates the need for any accessories.

All Oakmed pouches come in colour-coded packaging, which is easy to identify in storerooms and product guides.